

Texan™ Reagent Reservoirs for Multi-Well Pipettors

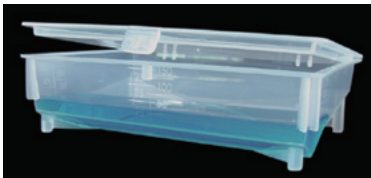


The design of our 175 mL Texan™ reagent reservoirs has been refined to include:

- volume graduations to 150 mL
- a reinforced base that won't wobble on the counter, avoiding spills
- a new easy-to-open latch that can be manipulated with one hand
- greater clarity for improved visibility of contents

Texan reservoirs are reusable – an environmentally sound alternative to disposable reservoirs. Their extra-large capacity means fewer refills and time saved. The V-shaped basin is compatible with both eight- and twelve-place pipettors with minimal reagent waste. Lids are removable – attached with ball and socket hinges – and latch tightly for secure storage of left-over reagents. With lids closed, Texan reservoirs are stackable. A perimeter ridge on the lid prevents sliding. Construction is clear, chemical-resistant, autoclavable polypropylene. Molded graduation marks provide volume indication at 25, 50, 75, 100, 125, and 150 mL.

Available sterile and non-sterile. Sterile reservoirs are supplied in bags of 5, nonsterile in bags of 10, both within outer packages of 100.

**Catalog No.****Description**

RSV-L175-S05

Texan Reservoirs with Lid, Pkg of 100, **Sterile**

RSV-0175-S05

Texan Reservoirs without Lid, Pkg of 100, **Sterile**

RSV-L175-N10

Texan Reservoirs with Lid, Pkg of 100, **Non-Sterile**

RSV-0175-N10

Texan Reservoirs without Lid, Pkg of 100, **Non-Sterile** **EXCEL Scientific, Inc.**
The Plate-Sealing and Reagent-Handling Experts


18350 George Blvd.
Victorville, CA 92394
email: info@excelscientific.com
760-246-4545 • FAX: 760-246-5588
www.excelscientific.com

"Committed to Exceeding Expectations!"

LINKwww.excelscientific.com/texan_content.html