

Zone-Free™


For easy sample access without adhesive fouling of pipet tips or probes

i Application:

Short-term storage and robotics

🌡️ Temperature Range:

-40°C to +90°C

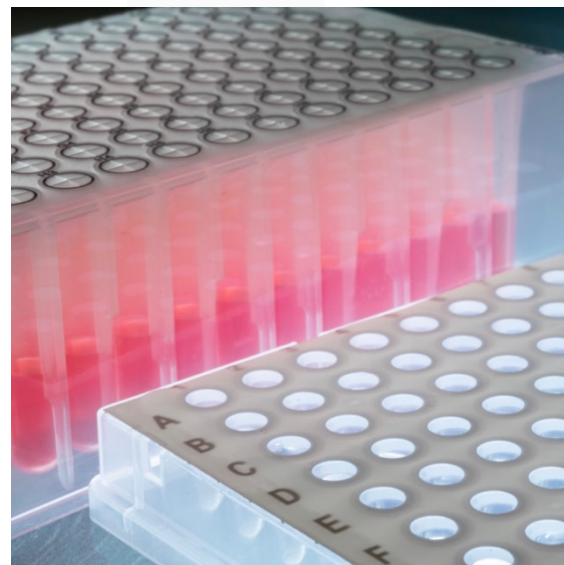

CATALOG #	DESCRIPTION	QUANTITY
ZAF-PE-50	Zone-Free	50 Sheets

Features and Benefits

- Offers excellent chemical resistance
- A multilayer film consisting of a top pierceable clear polyethylene film and a white polypropylene sublayer with 96 circular areas void of adhesive
- The white sublayer aids in the alignment of the film to the plate; alignment of opposing corner wells ensures that all 96 wells are properly aligned
- Adhesive-free zones are 9mm on center to facilitate multi-channel pipetting without adhesive fouling of pipet tips or robotic probes
- Pierceable
- Samples have minimal contact with adhesive due to the adhesive-free zones

Bundling Opportunities

- 96-well plates
- Pipettors
- Pipet tips and robotic probes
- Reservoirs